

Baumgartner, Adi ORCID:

<https://orcid.org/0000-0001-7042-0308>, Hartshorne, C., Polyzos, A., Weier, H.-U.G., Weier, J. F. and O'Brien, b. (2018) Full Karyotype Interphase Cell Analysis. Journal of Histochemistry and Cytochemistry, 66 (8). pp. 595-606.

Downloaded from: <http://ray.yorks.ac.uk/id/eprint/2932/>

The version presented here may differ from the published version or version of record. If you intend to cite from the work you are advised to consult the publisher's version:

Research at York St John (RaY) is an institutional repository. It supports the principles of open access by making the research outputs of the University available in digital form. Copyright of the items stored in RaY reside with the authors and/or other copyright owners. Users may access full text items free of charge, and may download a copy for private study or non-commercial research. For further reuse terms, see licence terms governing individual outputs. [Institutional Repository Policy Statement](#)

RaY

Research at the University of York St John

For more information please contact RaY at ray@yorks.ac.uk

Full Karyotype Interphase Cell Analysis

A. Baumgartner^{1, 2, 3 *}, C. Hartshorne¹, A. Polyzos², H.-U.G. Weier², J.F. Weier^{1,2,4}
and B. O'Brien^{2,5,6}

¹ Dept. of Obstetrics, Gynecology and Reproductive Sciences, University of California, San Francisco, CA, USA

² Life Sciences Division, E.O. Lawrence Berkeley National Laboratory, Berkeley, CA, USA

³ Biomedical Science, School of Health Sciences, York St John University, York, UK

⁴ Dermatopathology Service, University of California, San Francisco, CA, USA

⁵ Department of Perioperative Medicine, St Bartholomew's Hospital & Barts Heart Centre, London, UK

⁶ Outcomes Research Consortium, Cleveland Clinic, Cleveland, OH, USA

* To whom correspondence should be addressed:

Dr A Baumgartner, Biomedical Science, School of Health Sciences, York St John University,
York, YO31 7EX, UK

Phone: +44 190 487 6714, Email: a.baumgartner@yorks.ac.uk

Keywords: aneuploidy, full karyotype analysis, interphase nuclei, FISH, DNA probes, spectral imaging, chromosome enumeration

Running title: Full karyotype interphase cell analysis

Abstract

Aneuploidy seems to play not only a decisive role in embryonal development but also in tumorigenesis where chromosomal and genomic instability reflects a universal feature of malignant tumors. The cost of whole genome sequencing has fallen significantly, but it is still prohibitive for many institutions and clinical settings. No applied, cost-effective and efficient technique has been introduced yet aiming at research to assess the ploidy status of all 24 different human chromosomes in interphases simultaneously, especially in single cells. Here we present the selection of human probe DNA and a technique using multi-step FISH employing four sets of six labelled FISH probes able to delineate all 24 human chromosomes in interphase cells. This full karyotype analysis approach will provide additional diagnostic potential for single cell analysis. The use of spectral imaging (SI) has enabled the use of up to eight different fluorochrome labels simultaneously. Thus, scoring can be easily assessed by visual inspection, since SI permits computer-assigned and distinguishable pseudo-colors to each probe during image processing. This enables full karyotype analysis by FISH of single cell interphase nuclei.

Introduction

Evaluating numerical abnormalities and structural aberrations of chromosomes by fluorescence *in situ* hybridization (FISH) has been extensively carried out on cell metaphases and interphases since this technique has been introduced in the mid-1980s.¹ Methodologies such as spectral karyotyping (SKY)² or multiplex FISH³ and also comparative genomic hybridization (CGH)⁴ were able to effectively evaluate cytogenetic damage on metaphases across the whole genome. Nowadays, modern techniques like high-resolution chip-based CGH arrays and next-generation sequencing (NGS), in particular massive parallel sequencing (MPS), are capable of evaluating a plethora of cytogenetic changes. For array CGH, segmental DNA copy number variations at kilobase-pair resolution can be detected;⁵ while MPS is capable to analyze large parts of the genome by using shallow or low-pass whole genome sequencing when no coverage of the full genome is required, e.g. for preimplantation genetic diagnosis (PGD).⁶

In the early 2000s implantation rates had improved due to PGD aneuploidy screening by using commercially available chromosomal probe sets for single-cell analysis allowing the enumeration of up to ten chromosomes with fluorescence filter-based evaluation. Thereby, more than half of the numerical abnormalities seen in abnormal embryos originating from non-disjunction of chromosomes during cell division were covered.⁷⁻¹⁰ Without a doubt, aneuploidy is the most common cause of chromosomal abnormalities in humans leading to pregnancy loss.^{11, 12} Hence, shifting from partial karyotype FISH analysis to genome sequencing in the last decade allowed for simultaneous testing of numerous genetic aberrations and abnormalities. This was also evident when looking at the biopsied specimens, shifting away from polar bodies or blastomeres towards the trophectoderm.¹³ In human embryos, multiple molecular mechanisms that may also be involved in cancer formation can lead to aneuploidy and chromosomal mosaicism, thus to negative pregnancy outcomes; however, low-level mosaicism in human development may be a normal feature after all.¹⁴

Aneuploidy seems to also play a crucial role in cells of the extra-embryonal tissue that are important in implantation during early pregnancy and the formation of the placenta. These so-called invasive cytotrophoblasts (iCTBs) showed different aneuploidy levels on the basis

of their invasive behavior when assessed by using spectral imaging targeting six different chromosomes.¹⁵

Aneuploidy seems to play not only a decisive role in embryonal development but also in tumorigenesis where chromosomal and genomic instability reflects a universal feature of malignant tumors.¹⁶ It seems that the primary cause of pre-neoplastic / neoplastic genomic instability is the progression from stable diploid cells to unstable aneuploid cell species¹⁷ making aneuploidy a useful marker of malignant transformation.¹⁸

Although the cost for sequencing the whole genome has fallen to around \$1,000 per analyzed genome¹⁹ the price tag for equipment and material is quite high.²⁰ Even though larger hospital trusts, major universities and private biotech companies may have the funds to carry out high-throughput array chip methods and NGS on a daily basis, no applied, cost-effective and efficient technique has been introduced yet aiming at research to assess the ploidy status of all 24 different human chromosomes in interphases simultaneously, especially in single cells. Improving the coverage of all the chromosomes and devising sophisticated, fast and reliable methods to evaluate these single cells is favorable as not all of the possible numerical abnormalities can be currently assessed due to the limited number of available chromosome-specific probes and the limited number of suitable fluorochromes.^{21, 22}

Here we present the selection of human probe DNA and a technique using multi-step FISH employing four sets of six labelled FISH probes able to delineate all 24 human chromosomes in interphase cells. This full karyotype analysis approach will provide additional diagnostic potential for single cell analysis.

Material and Methods

Clone selection and DNA preparation

BAC clones²³ from the RP11 library (Invitrogen, Gaithersburg, MD) were chosen based on information available from the UC Santa Cruz (UCSC) genome sequence database (<http://genome.ucsc.edu/cgi-bin/hgGateway>) and the U.S. National Institute of Health, National Center for Biotechnology Information (NIH/NCBI) (http://www.ncbi.nlm.nih.gov/mapview/map_search.cgi?taxid=9606). Comprehensive DNA sequence information as well as structural organization of these BACs can be found in the above-mentioned databases. The preparation of DNA from BAC clones has been described in detail before.²⁴ Briefly, clones were cultured overnight in 10 ml Luria-Bertani (LB) medium containing 12.5 µg/ml chloramphenicol (Sigma, St. Louis, MO) and DNA was isolated using an alkaline lysis DNA extraction protocol.²⁵ In addition to BAC clones (see Table 1 for a complete overview) bacterial plasmid clones from the Weier lab at the LBL have been used. For the isolation of plasmids from clones RMC16L006²⁶, pBS444/7, pBS864, pBS1131, pBS8B/9, pBS239'-5', pBS609/51, 680TA-4 and W21R2-TA13 a commercial Kit from Qiagen was used on an overnight LB culture containing 30 µg/ml ampicillin (RMC and pBS clones) or kanamycin (TA clones), respectively. Except for the RMC16L006 plasmid DNA, all other plasmid DNA has been employed as templates in PCR reactions. The final probe sets can be seen in Table 2.

For PCR, 100 ng genomic DNA (Sigma), BAC or plasmid DNA was used as template for DNA amplification. PCR reactions (50 µl) were performed using 0.02 U/ml Taq Polymerase (Invitrogen) or JumpStart Taq polymerase (Sigma) in 1x PCR buffer (Invitrogen), 1.5 mM MgCl₂, 0.2 mM of each dNTP, and 0.6 mM of the forward and reverse primers (Qiagen; Alameda, CA). The chromosomes 1 and 6 alpha-satellite primers have been used as described previously.²⁷ The generation of chromosome 17 and 18 specific probes has also been previously published.²⁸ BlueScript primers WBS2 (ctc gga att aac cct cac taa agg) and WBS4 (gaa ttg taa tac gac tca cta tag) for the DNA amplification of alpha-satellite repeats were employed for chromosomes 8, 10, and 12. For chromosome 7 and 11, primers M13F / M13R²⁹ and WA8 (gat ggt agt agg ca[a/t] [c/g]t[c/a] aca gag) / WA9 (gat ggt agt agg cat c[a/c]c [a/c]aa g[a/t/c]a), respectively, have been used to amplify chromosome-specific DNA. For the amplification of chromosome 9 and 21 specific probe DNA a single primer was

in use for both, the satellite III primer W21R2 (caa acg tgc tca aag taa ggg aat g) and Jun15 (ccc aag ctt gca tgc gaa ttc), respectively. After an initial denaturation step of 2 min at 95 °C, 35 PCR cycles followed: denaturation at 95 °C for 40 sec, primer annealing at 54 °C for 1 min, and primer extension at 72 °C for 2.5 min. Ramp time was set to 30 sec for the first step followed by 1 min for the next two steps. A final step at 72 °C for 10 min concluded the PCR. PCR products were confirmed on a 2% agarose gel by applying 5 µl of the PCR reaction mixed with 1 µl of 0.4 g/ml sucrose solution.

DNA labeling and FISH

Using a commercially available kit (BioPrime Kit, Invitrogen, Gaithersburg, MD), random priming was employed to label all the BAC, plasmid and PCR-derived probe DNA. After initial testing employing indirect labels;³⁰⁻³² the random priming process was slightly modified to incorporate various fluorochrome dUTPs into the probe DNA: Cy5-dUTP and Cy5.5-dCTP (Amersham, Arlington Heights, IN), DEAC-dUTP (PerkinElmer, Waltham, MA) as well as SpectrumGreen-dUTP, SpectrumRed-dUTP and SpectrumOrange-dUTP (Vysis, Abbott Molecular Inc, Des Plaines, IL). Regarding the four sets of chromosome-specific probes, Table 2 gives also an overview of their direct fluorescent labels.

For FISH, labeled probe DNA is mixed with blocking DNA and concentrated via precipitation in ethanol. Salmon sperm DNA is added to block non-specific binding of the probe, and human COT1 DNA is added to block repetitive DNA sequences in the probes from binding to sites spread throughout several chromosomes/loci. The desired combination of labeled probe is mixed using 2-5 µl of each individual probe, depending on intensity of signal, as previously described.¹⁵

This is then combined with 1 µl of human Cot-1 DNA® (1 mg/ml, Invitrogen), 1 µl of salmon sperm DNA (10 mg/ml; Invitrogen), and 7 µl of the hybridization master mix (78.6% formamide, 14.3% dextran sulfate in 1.43x SSC, pH 7.0; 20x SSC is 3 M sodium chloride, 300 mM tri-sodium citrate) and thoroughly mixed and denatured at 76 °C for 10 min. The hybridization mixture was then pre-annealed by incubating at 37 °C for 30 min (allowing the Cot-1 DNA® to anneal to non-chromosome specific DNA repeats on the probes). In parallel, the metaphase slides prepared from phytohemagglutinin (PHA)-stimulated peripheral blood

lymphocytes from a karyotypically normal male³¹ were denatured for 3 min at 76 °C in 70% formamide/2x SSC, pH 7.0, dehydrated in 70%, 85%, and 100% ethanol for 2 min each, and allowed to air-dry. The hybridization mixture was then carefully applied to the slides, covered with a 22x22 mm² cover-slip and sealed with rubber cement. Slides were incubated overnight in a moist chamber at 37 °C. After removing rubber cement and the cover-slips, the slides were washed in 0.1x SSC at 43°C for 2 min, then, when biotin or digoxigenin labels have been used, incubated in PNM blocking reagent (5% nonfat dry milk powder (NESTLÉ Carnation, Wilkes-Barre, PA), 1% Nonidet-P40 (Sigma), 1% sodium azide (Sigma), 0.1 M sodium phosphate buffer, pH 8.0) for 10 min at room temperature. Bound probes were detected with fluorescein-conjugated avidin (avidin DCS, Vector labs, Burlingame, CA) and rhodamine-labeled anti-digoxigenin antibodies (Roche Diagnostics, Indianapolis, IN). In the case of direct-labeled probes, no immuno-detection step was necessary. Finally, after a last wash in PN or 2xSSC the slides were mounted with 4',6-diamidino-2-phenylindole (DAPI, 0.5 µg/ml; Calbiochem, La Jolla, CA) in antifade solution.^{32, 33}

Image acquisition and analysis

Fluorescence microscopy was performed on a Zeiss Axioskop microscope equipped with a SKY filter set (ChromaTechnology, Brattleboro, VT) for simultaneous observation of SKY suitable fluorochromes and also a DAPI filter (ChromaTechnology, Brattleboro, VT) for the detection of the counter stain. Images were collected using a cooled CCD camera (CCD-1300DS, VDS Vosskuehler, Osnabrück, Germany).²⁴ Further processing and printing of the images were done using the image processing software Adobe Photoshop (Adobe Systems Inc., San Jose, CA).

Results

Our probe sets have been constructed by choosing the most suitable probe DNA and fluorochrome so that each chromosome-specific probe within its set provides specificity and similar high efficiency. The probes used (see Tables 1 and 2 for a complete overview) are either locus-specific or repeat-specific (targeting alpha-satellites or satellite III). They are labeled with the following fluorescent labels: DEAC (excitation wavelength: 432 nm / emission wavelength: 472 nm, light blue), SpectrumGreen (497 nm / 524 nm, greens), SpectrumOrange (559 nm / 588 nm, orange), SpectrumRed (592 nm / 612 nm, red), Cy5 (649 nm / 666 & 670 nm, infrared), and Cy5.5 (675 nm / 694 nm, infrared).

Sets I and II detect chromosomes 13, 14, 16, 20, 21, and 22 as well as of chromosomes 15, 17, 18, 19, X, and Y, respectively (Figure 1) according to their fluorochrome labels (Figure 2), whereas sets III and IV are able to evaluate the rest of the chromosomes for the full karyotype by detecting chromosomes 2, 3, 4, 5, 9, and 12 as well as chromosomes 1, 6, 7, 8, 10, and 11, respectively (Figure 3). All the signals produced in metaphases and interphases by FISH are unambiguous, strong and do not cross-hybridize to other chromosomes (Figures 1 and 3). The exception was chromosome 9 in set III (Figure 3) which did produce a strong signal when tested for itself but was rather dim when co-hybridized with the rest of set III probes.

Discussion

Detection techniques for assessing numerical abnormalities and other cytogenetic aberrations often utilizes cost-effective and rapid classical staining methods such as Giemsa staining for karyotyping, but also molecular diagnostic tools like FISH and its multiple variations for a quick, robust and reliable detection of genetic damage. During the last decade, methodology development has been progressing towards a full karyotype analysis;³⁴ however, adjusting probe sets, i.e. probe target and fluorochrome label, rapidly to particular needs in the laboratory is still quite difficult. Using BAC clones as a source of probe DNA for FISH is cheap and effective^{28, 35, 36} on the other hand allows the analysis of the whole genome^{37, 38} for cytogenetic diagnoses. However, this metaphase-specific approach cannot be used for analysis of a single cells *per se*, as they are likely to be in an interphase stage.

While several bright chromosome-specific DNA repeat probes have been prepared and cloned by our labs in previous years,^{30, 39} the approach does not work for all human chromosomes. The acrocentric chromosomes 13 and 21, for example, share a high level homologous of sequences which is found heteromorphic in some individuals.⁴⁰ The BAC approach⁴¹ has advantages for complete chromosome enumeration in interphase cells. During preparation of this manuscript, Ioannou *et al.* demonstrated this using BAC clones from Roswell Park Cancer Institute RP-11 library. We decided to combine preexisting DNA repeat probes with optimized BAC contigs which identified using bioinformatics^{22, 42} to obtain optimal specificity and brightness (Table 2).

Now, combining the versatility of BACs and preexisting repeat probes with a wide repertoire of different fluorochromes together with the use of the SKY system resulted in a cost-effective, flexible and reliable methodology to detect four sets of six probes in interphase nuclei (see Figures 1 and 3). Incorporating consecutive hybridization steps (cell recycling) considerably increased the diagnostic range of existing FISH technologies.⁴³

Sequential multi-locus interphase FISH is one strategy where chromosome specific (locus-specific and alpha-centromeric) FISH probes have been hybridized sequentially to the same cell, initially done in formalin fixed and paraffin embedded histological sections from tumor tissues.⁴⁴ The use of spectral imaging (SI) has enabled the use of up to eight different

fluorochrome labels (with emission spectra from 450-1000 nm) simultaneously.⁴⁵ The scoring can be easily assessed by visual inspection, since SIm permits computer-assigned and distinguishable pseudo-colors to each probe during image processing.² This lead to the full karyotype analysis by FISH of single cell interphase nuclei like those of iCTBs in placental cell analysis (data not shown in this publication).

The benefits of combing sets of single copy DNA probes with separate sets of DNA repeat probes that contain super-bright signals in sequential hybridizations may raise concerns of DNA loss in repeated cycles of denaturation, hybridization and wash steps. The approach presented here does not eliminate potential problems associated with DNA losses, but optimized BAC contigs and plasmids targeting highly reiterated DNA repeat sequences consistently leads to brighter signals, thus reducing the negative effect of said losses.

As seen in Figure 3, the flexibility of choosing probe DNA and fluorochromes may have an unforeseen negative consequence, chromosome 9 showed a very dim signal which was below the threshold of detection for the SKY system, even though the same probe resulted in good strong signals when tested individually. Hence, either cell-type or donor differences as well as the interaction of DNA probes with each other within a set could lead to a less prominent of even a very dim and almost not visible signal. Thus, quality control when tailoring such probe sets to work in particular cell types is imperative.

In this publication, we presented the development of a probe collection made up of four sets of six labeled chromosome-specific FISH probes which can easily be modified towards three sets of eight probes using far-infrared fluorochromes like Cy7. These have been arranged so that the first two sets detect the most prevalent numerical abnormalities observed in human embryos, and the last two sets fill in the gap towards a full karyotype analysis. The generation of these probe sets shows the full potential of BAC/plasmid clones, which can be rapidly selected and tailored for specific genetic screening applications. The strong signal intensities from these repeat-rich probes and the labeling methodology employed, allows reduction in the costs for a single hybridization event by approximately 10-fold over conventional FISH.

In a clinical PGD setting, there is no apparent reason to do further analysis using multiple FISH probe sets after determining an abnormal ploidy status such as trisomy 13; however, if information on all 24 different chromosomes for individual embryos, i.e. the discarded chromosomally abnormal embryos, can be collected, this could be very important information for early embryo development, and this information may be useful for future clinical diagnosis. Also, with regard to array and NGS analysis in a clinical PGD setting, 250-500 ng of DNA, an equivalent of about 35,000-70,000 cells, would be required for analysis. This entails that if there is only a limited number of cells (i.e. blastocytes) whole genome amplification will have to be applied increasing time and cost but also errors.⁴⁶ As our 24-probe set was originally developed for and applied to score all 24 chromosomes in placental invasive cytotrophoblasts (CTBs), we have found that at least 75% of the male invasive CTBs have gained extra copies of chromosomes with the most common gains being acrocentric chromosomes. Also, most of the invasive CTBs did not have the same karyotype (unpublished data). The probes that have been developed and tested are typically very useful to do full karyotyping on the few interphase cells, such as cancer stem cells, fetal cells in maternal blood, or heterogeneous cells with different karyotypes.

The set of 4x6 chromosome enumerating FISH probes has been initially developed in-house to study invasive placental cells and to aid in PGD due to the fact that only some commercial probes are available as enumeration sets with a maximum of four different colors. There are many clinical applications for DNA probe sets like the ones described here, i.e. prenatal analysis of aneuploidy. But for example, if one wishes to study malignant mesothelioma, the panel would include a probe for *ERBB2*.⁴⁷ Similarly, overexpression due to potential amplifications of the epidermal growth factor receptor (*EGFR*) gene in squamous carcinoma cell lines might include an EGF receptor-specific DNA probe.⁴⁸ Hence, it is very conceivable that our set of probes, original or modified, could find a use in other studies, such as in tumor cytogenetic evaluations or chromosome analysis in genotoxicology and mutagenesis, or whenever quick chromosome analysis of numerical abnormalities in interphase cells is required.

Acknowledgements & Funding Statement

This work was supported in part by NIH grants CA-88258, HD-045736 and ND-044313, and a grant from the Director, Office of Energy Research, Office of Health and Environmental Research, U.S. Department of Energy, under contract DE-AC02-05CH11231. A.B. was supported in part by a grant from the University of California Discovery Program (BIO03-10414 to JF). J.F. was supported by NIH grant HD-041425.

Disclaimer

This document was prepared as an account of work sponsored by the United States Government. While this document is believed to contain correct information, neither the United States Government nor any agency thereof, nor The Regents of the University of California, nor any of their employees, makes any warranty, express or implied, or assumes any legal responsibility for the accuracy, completeness, or usefulness of any information, apparatus, product, or process disclosed, or represents that its use would not infringe privately owned rights. Reference herein to any specific commercial product, process, or service by its trade name, trademark, manufacturer, or otherwise, does not necessarily constitute or imply its endorsement, recommendation, or favoring by the United States Government or any agency thereof, or The Regents of the University of California. The views and opinions of authors expressed herein do not necessarily state or reflect those of the United States Government or any agency thereof or The Regents of the University of California.

Competing Interest Statement

The authors declare they have no competing interests.

Author Contribution Statement

UW & BO designed the study, AB & UW created and optimized the DNA probes using BACs, plasmids and PCR; AB, CH, AP & JW optimized fluorescent labels and carried out fluorescence in situ hybridization and spectral imaging as well as analysis of single probes and set of probes, AB, BO, UW & JW drafted the manuscript. All authors have read and approved the final manuscript.

References

1. Pinkel D, Straume T and Gray JW. Cytogenetic analysis using quantitative, high-sensitivity, fluorescence hybridization. *Proc Natl Acad Sci U S A* 1986; 83: 2934-2938.
2. Schrock E, du Manoir S, Veldman T, Schoell B, Wienberg J, Ferguson-Smith MA, Ning Y, Ledbetter DH, Bar-Am I, Soenksen D, Garini Y and Ried T. Multicolor spectral karyotyping of human chromosomes. *Science* 1996; 273: 494-497.
3. Speicher MR, Gwyn Ballard S and Ward DC. Karyotyping human chromosomes by combinatorial multi-fluor FISH. *Nat Genet* 1996; 12: 368-375. DOI: 10.1038/ng0496-368.
4. Kallioniemi A, Kallioniemi OP, Sudar D, Rutovitz D, Gray JW, Waldman F and Pinkel D. Comparative genomic hybridization for molecular cytogenetic analysis of solid tumors. *Science* 1992; 258: 818-821.
5. Le Scouarnec S and Gribble SM. Characterising chromosome rearrangements: recent technical advances in molecular cytogenetics. *Heredity (Edinb)* 2012; 108: 75-85. DOI: 10.1038/hdy.2011.100.
6. Deleye L, Dheedene A, De Coninck D, Sante T, Christodoulou C, Heindryckx B, Van den Abbeel E, De Sutter P, Deforce D, Menten B and Van Nieuwerburgh F. Shallow whole genome sequencing is well suited for the detection of chromosomal aberrations in human blastocysts. *Fertil Steril* 2015; 104: 1276-1285 e1271. DOI: 10.1016/j.fertnstert.2015.07.1144.
7. Baart EB, Martini E and Van Opstal D. Screening for aneuploidies of ten different chromosomes in two rounds of FISH: a short and reliable protocol. *Prenat Diagn* 2004; 24: 955-961. DOI: 10.1002/pd.1052.
8. Baart EB, Van Opstal D, Los FJ, Fauser BC and Martini E. Fluorescence in situ hybridization analysis of two blastomeres from day 3 frozen-thawed embryos followed by analysis of the remaining embryo on day 5. *Hum Reprod* 2004; 19: 685-693. DOI: 10.1093/humrep/deh094.
9. Bahce M, Escudero T, Sandalinas M, Morrison L, Legator M and Munne S. Improvements of preimplantation diagnosis of aneuploidy by using microwave hybridization, cell recycling and monocolour labelling of probes. *Mol Hum Reprod* 2000; 6: 849-854.
10. Munne S, Magli C, Bahce M, Fung J, Legator M, Morrison L, Cohert J and Gianaroli L. Preimplantation diagnosis of the aneuploidies most commonly found in spontaneous abortions and live births: XY, 13, 14, 15, 16, 18, 21, 22. *Prenat Diagn* 1998; 18: 1459-1466.
11. Hassold T, Hall H and Hunt P. The origin of human aneuploidy: where we have been, where we are going. *Hum Mol Genet* 2007; 16 Spec No. 2: R203-208. DOI: 10.1093/hmg/ddm243.

12. Simpson JL. Causes of fetal wastage. *Clin Obstet Gynecol* 2007; 50: 10-30. DOI: 10.1097/GRF.0b013e31802f11f6.
13. Treff NR and Zimmerman RS. Advances in Preimplantation Genetic Testing for Monogenic Disease and Aneuploidy. *Annu Rev Genomics Hum Genet* 2017; 18: 189-200. DOI: 10.1146/annurev-genom-091416-035508.
14. McCoy RC. Mosaicism in Preimplantation Human Embryos: When Chromosomal Abnormalities Are the Norm. *Trends Genet* 2017; 33: 448-463. DOI: 10.1016/j.tig.2017.04.001.
15. Weier JF, Hartshorne C, Nguyen HN, Baumgartner A, Polyzos AA, Lemke KH, Zeng H and Weier HU. Analysis of human invasive cytotrophoblasts using multicolor fluorescence in situ hybridization. *Methods* 2013; 64: 160-168. DOI: 10.1016/j.ymeth.2013.05.021.
16. Todorovic-Rakovic N. Genome-based versus gene-based theory of cancer: possible implications for clinical practice. *J Biosci* 2011; 36: 719-724.
17. Duesberg P, Fabarius A and Hehlmann R. Aneuploidy, the primary cause of the multilateral genomic instability of neoplastic and preneoplastic cells. *IUBMB Life* 2004; 56: 65-81. DOI: 10.1080/15216540410001667902.
18. Alaizari NA, Sperandio M, Odell EW, Peruzzo D and Al-Maweri SA. Meta-analysis of the predictive value of DNA aneuploidy in malignant transformation of oral potentially malignant disorders. *J Oral Pathol Med* 2017. DOI: 10.1111/jop.12603.
19. Hayden EC. Technology: The \$1,000 genome. *Nature* 2014; 507: 294-295. DOI: 10.1038/507294a.
20. Liehr T. What about the real costs of next generation sequencing (NGS) in human genetic diagnostics?, <https://atlasofscience.org/what-about-the-real-costs-of-next-generation-sequencing-ngs-in-human-genetic-diagnostics/> (2017, accessed March 2018).
21. Lemke KH, Weier JF, Weier HG and Lawin-O'Brien AR. High Performance DNA Probes for Perinatal Detection of Numerical Chromosome Aberrations. *Adv Tech Biol Med* 2015; 3. DOI: 10.4172/2379-1764.1000155.
22. Hsu JH, Zeng H, Lemke KH, Polyzos AA, Weier JF, Wang M, Lawin-O'Brien AR, Weier HU and O'Brien B. Chromosome-specific DNA repeats: rapid identification in silico and validation using fluorescence in situ hybridization. *Int J Mol Sci* 2012; 14: 57-71. DOI: 10.3390/ijms14010057.
23. Shizuya H, Birren B, Kim UJ, Mancino V, Slepak T, Tachiiri Y and Simon M. Cloning and stable maintenance of 300-kilobase-pair fragments of human DNA in Escherichia coli using an F-factor-based vector. *Proc Natl Acad Sci U S A* 1992; 89: 8794-8797.

24. Weier HU, Tuton TB, Ito Y, Chu LW, Lu CM, Baumgartner A, Zitzelsberger HF and Weier JF. Molecular cytogenetic characterization of chromosome 9-derived material in a human thyroid cancer cell line. *Cytogenet Genome Res* 2006; 114: 284-291. DOI: 10.1159/000094215.
25. Birnboim HC and Doly J. A rapid alkaline extraction procedure for screening recombinant plasmid DNA. *Nucleic Acids Res* 1979; 7: 1513-1523.
26. Moyzis RK, Albright KL, Bartholdi MF, Cram LS, Deaven LL, Hildebrand CE, Joste NE, Longmire JL, Meyne J and Schwarzacher-Robinson T. Human chromosome-specific repetitive DNA sequences: novel markers for genetic analysis. *Chromosoma* 1987; 95: 375-386.
27. Dunham I, Lengauer C, Cremer T and Featherstone T. Rapid generation of chromosome-specific aliphoid DNA probes using the polymerase chain reaction. *Hum Genet* 1992; 88: 457-462.
28. Baumgartner A, Weier JF and Weier HU. Chromosome-specific DNA repeat probes. *J Histochem Cytochem* 2006; 54: 1363-1370. DOI: 10.1369/jhc.6A6974.2006.
29. O'Shaughnessy JB, Chan M, Clark K and Ivanetich KM. Primer design for automated DNA sequencing in a core facility. *Biotechniques* 2003; 35: 112-116, 118-121.
30. Fung J, Hyun W, Dandekar P, Pedersen RA and Weier HU. Spectral imaging in preconception/preimplantation genetic diagnosis of aneuploidy: multicolor, multichromosome screening of single cells. *J Assist Reprod Genet* 1998; 15: 323-330.
31. Fung J, Weier HU and Pedersen RA. Detection of structural and numerical chromosome abnormalities in interphase cells using spectral imaging. *J Histochem Cytochem* 2001; 49: 797-798. DOI: 10.1177/002215540104900616.
32. Weier HU, Wang M, Mullikin JC, Zhu Y, Cheng JF, Greulich KM, Bensimon A and Gray JW. Quantitative DNA fiber mapping. *Hum Mol Genet* 1995; 4: 1903-1910.
33. Weier HU, Rhein AP, Shadravan F, Collins C and Polikoff D. Rapid physical mapping of the human trk protooncogene (NTRK1) to human chromosome 1q21-q22 by P1 clone selection, fluorescence in situ hybridization (FISH), and computer-assisted microscopy. *Genomics* 1995; 26: 390-393.
34. Turner K, Fowler K, Fonseka G, Griffin D and Ioannou D. Multicolor detection of every chromosome as a means of detecting mosaicism and nuclear organization in human embryonic nuclei. *Panminerva Med* 2016; 58: 175-190.
35. Lu CM, Kwan J, Baumgartner A, Weier JF, Wang M, Escudero T, Munne S, Zitzelsberger HF and Weier HU. DNA probe pooling for rapid delineation of chromosomal breakpoints. *J Histochem Cytochem* 2009; 57: 587-597. DOI: 10.1369/jhc.2009.953638.

36. Weier HU, Ito Y, Kwan J, Smida J, Weier JF, Hieber L, Lu CM, Lehmann L, Wang M, Kassabian HJ, Zeng H and O'Brien B. Delineating chromosomal breakpoints in radiation-induced papillary thyroid cancer. *Genes (Basel)* 2011; 2: 397-419. DOI: 10.3390/genes2030397.
37. Voullaire L, Slater H, Williamson R and Wilton L. Chromosome analysis of blastomeres from human embryos by using comparative genomic hybridization. *Hum Genet* 2000; 106: 210-217.
38. Wells D and Delhanty JD. Comprehensive chromosomal analysis of human preimplantation embryos using whole genome amplification and single cell comparative genomic hybridization. *Mol Hum Reprod* 2000; 6: 1055-1062.
39. Weier HU, Kleine HD and Gray JW. Labeling of the centromeric region on human chromosome 8 by in situ hybridization. *Hum Genet* 1991; 87: 489-494.
40. Weier HU and Gray JW. A degenerate alpha satellite probe, detecting a centromeric deletion on chromosome 21 in an apparently normal human male, shows limitations of the use of satellite DNA probes for interphase ploidy analysis. *Anal Cell Pathol* 1992; 4: 81-86.
41. Ioannou D, Meershoek EJ, Thornhill AR, Ellis M and Griffin DK. Multicolour interphase cytogenetics: 24 chromosome probes, 6 colours, 4 layers. *Mol Cell Probes* 2011; 25: 199-205. DOI: 10.1016/j.mcp.2011.08.002.
42. Zeng H, Weier JF, Wang M, Kassabian HJ, Polyzos AA, Baumgartner A, O'Brien B and Weier HU. Bioinformatic Tools Identify Chromosome-Specific DNA Probes and Facilitate Risk Assessment by Detecting Aneusomies in Extra-embryonic Tissues. *Curr Genomics* 2012; 13: 438-445. DOI: 10.2174/138920212802510510.
43. Harrison CJ, Gibbons B, Yang F, Butler T, Cheung KL, Kearney L, Dirscherl L, Bray-Ward P, Gregson M and Ferguson-Smith M. Multiplex fluorescence in situ hybridization and cross species color banding of a case of chronic myeloid leukemia in blastic crisis with a complex Philadelphia translocation. *Cancer Genet Cytogenet* 2000; 116: 105-110.
44. Walch A, Bink K, Hutzler P, Bowering K, Letsiou I, Zitzelsberger H, Braselmann H, Stein H, Hofler H and Werner M. Sequential multilocus fluorescence in situ hybridization can detect complex patterns of increased gene dosage at the single cell level in tissue sections. *Lab Invest* 2001; 81: 1457-1459.
45. Weier JF, Ferlatte C, Baumgartner A, Jung CJ, Nguyen HN, Chu LW, Pedersen RA, Fisher SJ and Weier HU. Molecular cytogenetic studies towards the full karyotype analysis of human blastocysts and cytotrophoblasts. *Cytogenet Genome Res* 2006; 114: 302-311. DOI: 10.1159/000094218.

46. Sho S, Court CM, Winograd P, Lee S, Hou S, Graeber TG, Tseng HR and Tomlinson JS. Precision oncology using a limited number of cells: optimization of whole genome amplification products for sequencing applications. *BMC Cancer* 2017; 17: 457. DOI: 10.1186/s12885-017-3447-6.
47. Horvai AE, Li L, Xu Z, Kramer MJ, Jablons DM and Treseler PA. Malignant mesothelioma does not demonstrate overexpression or gene amplification despite cytoplasmic immunohistochemical staining for c-Erb-B2. *Arch Pathol Lab Med* 2003; 127: 465-469. DOI: 10.1043/0003-9985(2003)127<0465:MMDNDO>2.0.CO;2.
48. Yamamoto T, Kamata N, Kawano H, Shimizu S, Kuroki T, Toyoshima K, Rikimaru K, Nomura N, Ishizaki R, Pastan I, Gamou S and Shimizu N. High incidence of amplification of the epidermal growth factor receptor gene in human squamous carcinoma cell lines. *Cancer Res* 1986; 46: 414-416.

Table legends

Table 1 legend

Table 1 shows a complete overview of BAC clones and bacterial plasmid clones used.

Table 2 legend

Table 2 shows the final probe set used.

It also gives an overview of the direct fluorescent labels for the four sets of chromosome-specific probes.

Table 1

Locus	Probe name	Clones/DNA	End sequence 1	End sequence 2	Full sequence ID
13q21.31	OR7E156P	RP11-527N12	-	-	AL354810
		RP11-282D7	-	-	AL355609
		RP11-320N6	-	-	AL359208
		RP11-67L17	-	-	AL354774
		RP11-473M10	-	-	AL445989
		RP11-394A14	-	-	AL445238
		RP11-520F9	-	-	AL355879
		RP11-205B18	-	-	AL354736
14q13.3	PAX9	RP11-12H15	B75808	B75809	-
		RP11-150O18	AQ378665	AQ378667	-
		RP11-452H6	AQ583099	AQ583102	-
		RP11-381L10	AQ532441	AQ554647	-
		RP11-73H19	AQ266602	AQ266604	-
		RP11-49P15	AQ051953	AQ051955	-
		RP11-151J2	AQ379285	AQ379286	-
		RP11-410J4	AQ549717	-	-
16qh, sat II	RMC16L006	RMC16L006	-	-	X06138
20p11.1-11.2	SRCext	RP11-298O1	AQ507400	AQ507403	-
		RP11-465M13	AQ636482	-	-
		RP11-192N1	AQ412321	AQ412322	-
		RP11-151C5	AQ376308	AQ376305	-
		RP11-451G10	AQ583252	AQ583256	-
		RP11-103B19	AQ313159	AQ313162	-
		RP11-467A7	AQ637700	AZ516714	-
		RP11-99B19	AQ318386	AQ318387	-
21q22	D21S167	YAC 141G6	-	-	X52289
22q11.22-q11.23	BCR	RP11-357H16	AZ518881	AQ553050	-
		RP11-62K15	AQ199674	AQ199676	-
		RP11-164N13	AQ380113	AQ380117	-
		RP11-116G21	AQ348695	AQ348692	-
		RP11-113C11	AQ344985	AQ344986	-
15q25.3	NTRK3ext	RP11-96B23	AQ313684	AQ313681	-
		RP11-114I9	AQ344858	AQ344856	-
		RP11-285I14	-	-	AC011966
		RP11-427O16	-	-	AC023844
		RP11-356B18	-	-	AC009711
		RP11-247E14	-	-	AC087593
		RP11-97O12	-	-	AC013489
		RP11-285M22	-	-	AC131274
		genomic DNA	-	-	M65181
17cen, α -sat	17cen	CTD-2052L21	AQ270406	AQ235108	-
		CTD-2288H11	B98832	-	-
		CTD-2195C23	-	-	-
		CTD-2017C04	-	-	-
		CTD-2195B23	-	-	AC011510
18cen, α -sat	18cen	CTD-2218D21	AQ032786	-	-
		RP11-422C9	AQ553389	AQ553391	-
		RP11-475A12	AQ635998	AQ636000	-
		RP11-14A18	B76225	B82244	-
		RP11-242E13	-	-	AC068123
19q13.2	AXL	genomic DNA	-	-	-
Xq21	CYCL1	RP11-65A9	AQ237129	AQ237131	-
		RP11-139C21	AQ382574	AQ382578	-
		RP11-71D7	AQ236987	AQ267891	-
		RP11-373L24	-	-	AC010733
Yq12	RP11-242E13	genomic DNA	-	-	-
1cen, α -sat	1cen				
2p16.1-15	RELex				

		RP11-375M18	AQ533441	AQ551245	-
		RP11-418N22	AQ550069	AQ550072	-
		RP11-77P21	AQ284573	AQ284575	-
		RP11-477N2	AQ637330	AQ637326	-
3q27.3	BCL6	RP11-143D11	-	-	AC092103
		RP11-567G11	-	-	AC104635
		RP11-88P6	-	-	AC018919
		RP11-211G3	-	-	AC072022
		RP11-58M14	AQ199229	AQ199231	-
		RP11-120O8	AQ350519	AQ350515	-
4q22.1	TIGD2ext	RP11-1E24	B48349	AQ312932	-
		RP11-502A23	-	-	AC079141
		RP11-84C13	-	-	AC104785
		RP11-173C9	-	-	AC105388
		RP11-549C16	-	-	AC093862
		RP11-15F14	B76416	B76417	-
		RP11-115D19	-	-	AC097478
		RP11-67M1	-	-	AC093759
		RP11-350B19	-	-	AC105445
		RP11-176N15	-	-	AC108038
5q23.1	05BP1-S47	RP11-183D16	-	-	AC093781
		RP11-23E11	B86433	B86434	-
		RP11-254M1	AQ479016	AQ479018	-
		RP11-464H3	AQ586366	AZ515952	-
		RP11-133L2	AQ350910	AQ350911	-
		RP11-59G17	AQ194868	AQ194864	-
		RP11-42O22	AQ116158	AQ046673	-
		RP11-185N19	-	-	AC021224
12cen, α -sat	12cen	444/7	-	-	G03348
6cen, α -sat	6cen	864	-	-	G04505
7cen, α -sat	7cen	680TA-4	-	-	AJ295152
8cen, α -sat	8cen	8B/9	-	-	M64779
9cen, sat III	9cen	W21R2-TA13	-	-	X06137
10cen, α -sat	10cen	609/51	-	-	X63622
11cen, α -sat	11cen	238'-5'	-	-	M21452

Table 2

Set	Locus	Position (Mbp)	Probe type	Label	Color	Clone/Contig ID
I	13q21.31	62.520-63.638	BAC pool of 8 BACs	Cy5	IR1	OR7E156P
I	14q13.3	35.678-36.975	BAC pool of 8 BACs	DEAC	aqua	PAX9
I	16qh, sat II	-	Plasmid	Cy5.5	IR2	RMC16L006
I	20q11.1-11.2	34.950-36.152	BAC pool of 9 BACs	Sp.Red	red	SRCext
I	21q22	around 37.770	PCR product	Sp.Green	green	D21S167
I	22q11.22/23	21.545-22.085	BAC pool of 3 BACs	Sp.Orange	orange	BCR
II	15q25.3-26.1	85.745-86.975	BAC pool of 7 BACs	Sp.Orange	orange	NTRK3ext
II	17cen, α -sat	-	PCR product	Sp.Red	red	17cen
II	18cen, α -sat	-	PCR product	Cy5.5	IR2	18cen
II	19q13.2	around 46.500	BAC pool of 6 BACs	Sp.Green	green	AXL
II	Xq21	82.447-82.915	BAC pool of 3 BACs	Cy5	IR1	CYCL1
II	Yq12	57.158-57.256	Single BAC	DEAC	aqua	RP11-242E13
III	2p16.1-15	60.525-61.831	BAC pool of 9 BACs	Sp.Red	red	RELext
III	3q27.3	188.590-188.976	BAC pool of 6 BACs	Sp.Green	green	BCL6
III	4q22.1	90.168-91.675	BAC pool of 10 BACs	Sp.Orange	orange	TIGD2ext
III	5q23.1	116.279-117.541	BAC pool of 7 BACs	DEAC	aqua	05BP1-S47
III	9cen, sat III	-	PCR product	Cy5.5	IR2	9cen
III	12cen, α -sat	-	PCR product	Cy5	IR1	12cen
IV	1cen, α -sat	-	PCR product	Cy5.5	IR2	1cen
IV	6cen, α -sat	-	PCR product	Sp.Green	green	6cen
IV	7cen, α -sat	-	PCR product	Sp.Red	red	7cen
IV	8cen, α -sat	-	PCR product	DEAC	aqua	8cen
IV	10cen, α -sat	-	PCR product	Sp.Orange	red	10cen
IV	11cen, α -sat	-	PCR product	Cy5	IR1	11cen

Figure legends

Figure 1

(A) FISH probe set I (see Table 2) hybridized on one metaphase spread from normal male lymphocyte and (B) one interphase nucleus. Computer assigned pseudo-colors can be seen showing two copies each of chromosomes 13, 14, 16, 20, 21, and 22. (C) FISH probe set II (see Table 2) hybridized on one metaphase spread from normal male lymphocyte and (D) one interphase nucleus. It showed two copies each of chromosomes 15, 17, 18, 19, and one copy of chromosome X and Y. The size bars on panels A and C indicate 5 μm while those on panels B and D indicate 2.5 μm .

Figure 2

(A) The emission spectra of DEAC, Spectrum Green, Spectrum Orange, Spectrum Red, Cy5, and Cy5.5 can be seen in this graph. (B) By using the distinguished emission spectra of these fluorochromes saved in a classified file the SKY system can easily identify individual chromosomes and create a karyotype. This is an example of an abnormal female cell hybridized with probe set II showing four copies each of chromosomes 15 and 19, three copies of chromosome X, and two copies each of chromosomes 17 and 18.

Figure 3

(A) FISH probe set III (see Table 2) hybridized on one normal male interphase nucleus (RGB colors). (B) The corresponding classified image from SKY system (pseudo-colors) showed two copies each of chromosomes 2, 3, 4, 5, and 12. The chromosome 9-Cy5.5 probe developed in-house showed weak hybridization signals and was not detected by the SKY system. (C) FISH probe set IV (see Table 2) hybridized on one normal male interphase nucleus (RGB colors). (D) The corresponding classified image from SKY system (pseudo-colors) showed two copies each of chromosomes 1, 6, 7, 8, 10, and 11. The size bar on panel A indicates 2.5 μm and is representative of all panels.

Figures

Figure 1

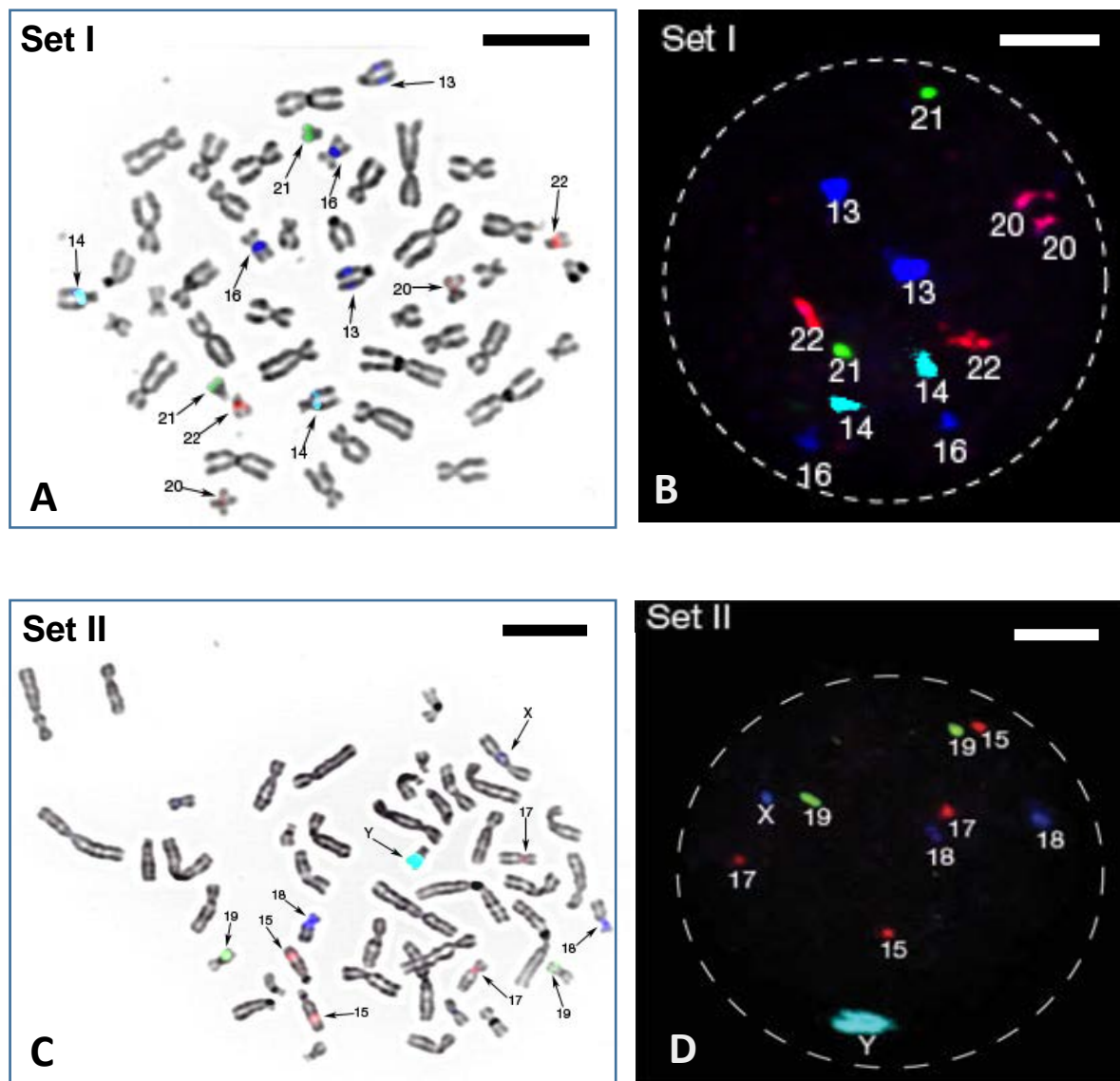


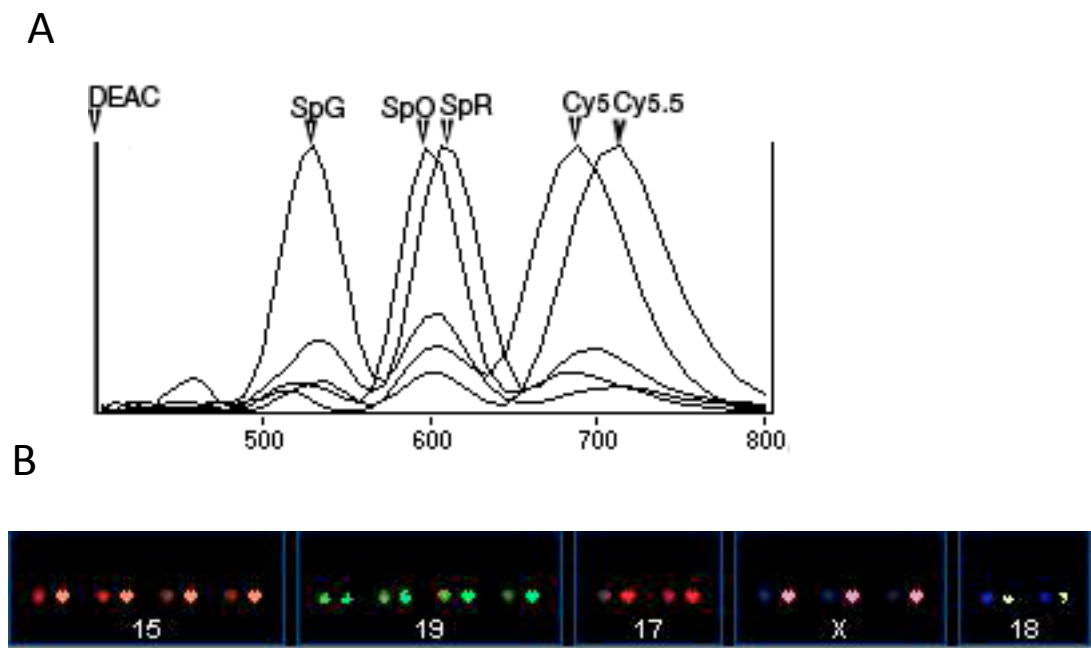
Figure 2

Figure 3

